




ZirChrom®



Technical
Bulletin #106

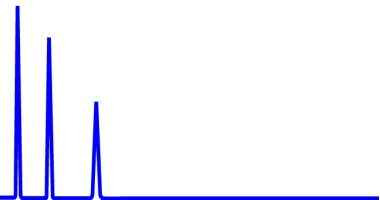

... For Peak Performance

ZirChrom®-PBD Stability

ZirChrom®-PBD columns are extremely stable over extended periods of use. They may be used at elevated column temperatures and over the entire pH range.

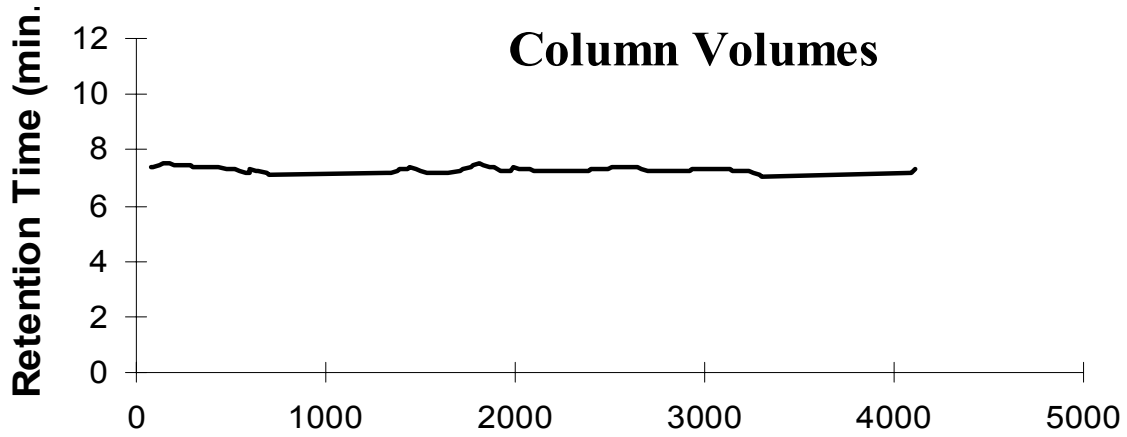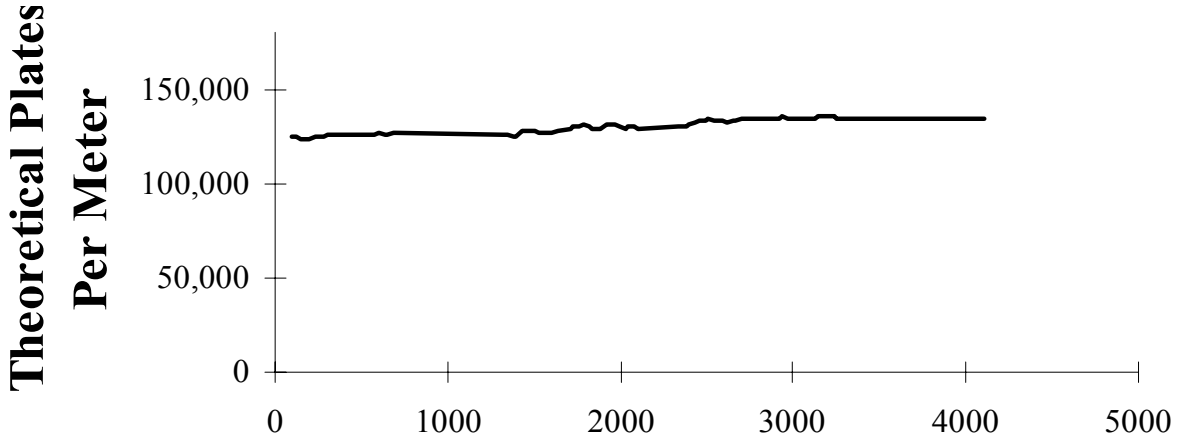

Column Volumes

LC Conditions

Column: ZirChrom®-PBD 150 × 4.6 mm
 Mobile Phase: 40/60 A/B
 A: Acetonitrile
 B: Water

Flow rate: 1.0 ml/min.
 Solute: Ethylbenzene

