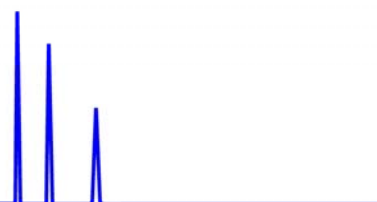




ZirChrom®

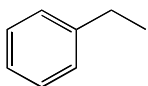
Technical
Bulletin #134

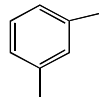


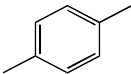
... For Peak Performance

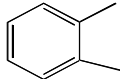
Isomer Separation on ZirChrom®-CARB

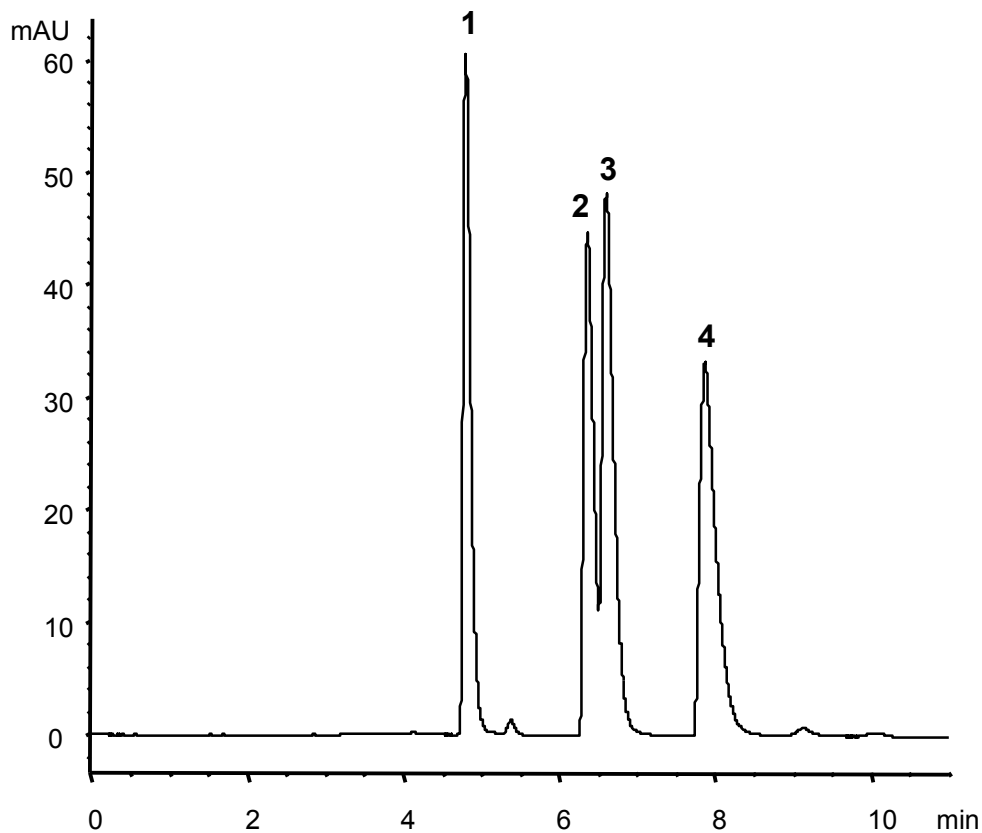
Analytes

1 - Ethylbenzene


2 - m-Xylene


3 - p-Xylene


4 - O-Xylene




LC Conditions

Column: A: ZirChrom®-CARB, 150 × 4.6 mm;

Mobile Phase: 50/50 A/B

A: Acetonitrile

B: Water

Flow rate: 1.0 mL/min.

Temperature: 30 °C

Detection: 254 nm

